

T304 1.5"—4" Dual Purge Port Ball Valves





Sanitary Solutions stocks a full port, cavity filled, T304 clamp connection dual purge port ball valve. Our dual purge port ball valve is constructed with a full port positioned inside the valve body between two polished flanges and two Teflon® valve seats. Take a moment to review the two purge ports located at the bottom of the valve for full drain ability, taking samples, and CIP/SIP applications.


Our dual purge port ball valve is 100% full opening, and equipped to take a maximum pressure of 200 PSI and temperature from 60°F to 300°F as well as maintaining a full flow rate with a minimum pressure drop.

In case of a tight situation or if your valve is close to a wall and will not open or close, you can remove the handle and re-position to open or close the valve. Please go to the www.sanitarysolutions.com for more information or contact us at 1-888-909-3569 or 1-855-640-7871.

www.sanitarysolutionsinc.com



T304 Dual Purge Port Ball Valve Dimensions

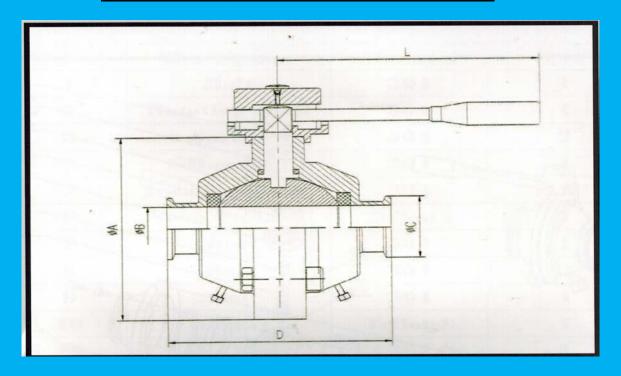

Tube O.D Size	<u>D</u>	<u>A</u>	<u>B</u>	<u>C</u>	L
1.5"	3.937"	4.528"	1.370"	1.988"	6.299"
2.0"	4.331"	5.157"	1.870"	2.520"	7.480"
2.5"	5.118"	6.299"	2.370"	3.051"	8.268"
3.0"	5.512"	7.323"	2.870"	3.583"	9.843"
4.0"	6.692"	8.661"	3.835"	4.685"	12.60"

"DUAL Purge Port Ball Valve" PART NUMBER EXAMPLE WHEN ORDERING: SSPPBV-SIZE-304

1-888-909-3569 or 1-855-640-7871



T304 Dual Purge Port Ball Valve Product Breakdown

Item	Description	Material	Qty.
S 2	Bolt	T304	1 1
2	Rotating Shaft	T304	1
3	Spring	T304	2
	Handle	T304	1
5	Locating Sleeve	T304	1
6	Stem Packing	Teflon®	1
7	Stem	T304	1
8	Bolt	T304	3
9	Nut	T304	
10	Washer	Teflon®	1
11	Valve Body	T304	1
12	O-Ring Seal	Buna	1
13	Seat	Teflon®	1
14	Gasket	T304	1
15	Ball Seat	T304	1
16	Seat	Teflon®	1
17	O-Ring Seal	Buna	1
18	Valve Body	T304	1
19	Bolt	T304	3
20	Washer	Teflon®	1
21	Nut	T304	1
22	Hex Deck Bolt	T304	1
23	Hex Deck Bolt	T304	1

